
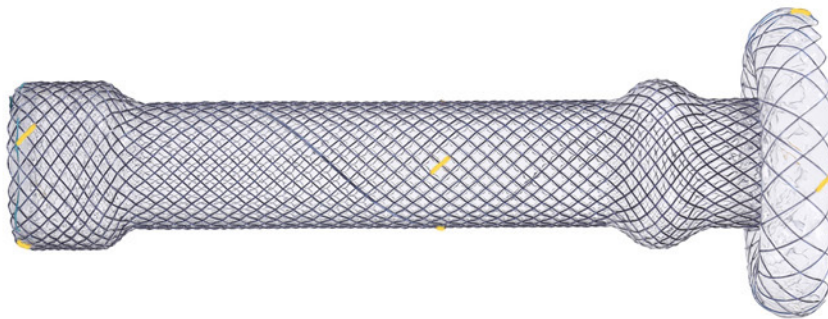


# Cardia Umbrella Fistula Occluding Esophageal Stent


**Application:** Esophageal or gastric anastomotic fistula

## Performance Characteristics:

- Umbrella for fistula-occluding
- Middle ball for better migration resistance
- Elastic covering to resist tumor inner growth
- Spherical part prevent side leakage effectively
- Distal and proximal deployment
- Delivery system with markers for better position accuracy

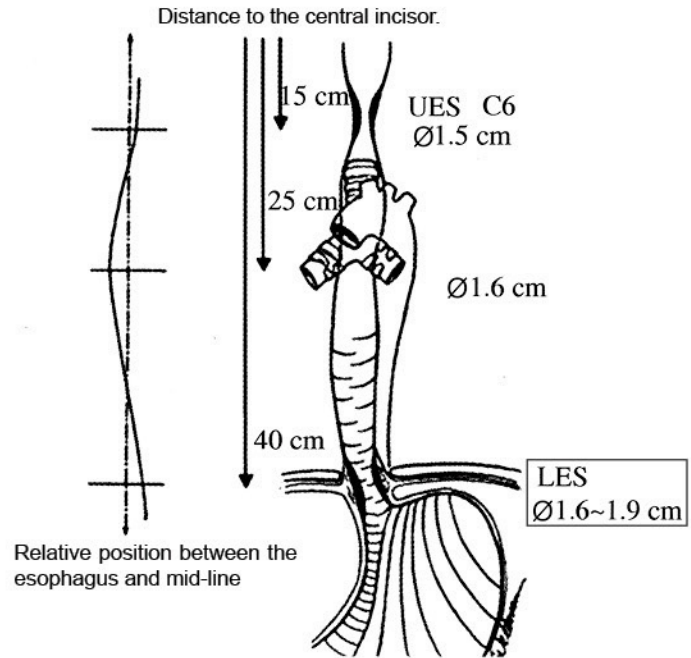

## Specifications:

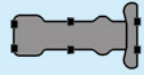Product Code	Dia. A (mm)	Total Length B (mm)	Product Illustration	Introducer
NTJB2-AB- HC/FC	A:18, cup and ball:24, umbrella:45	80/100/120/140		24Fr*700mm OTW
	A:20, cup and ball:26, umbrella:50			18Fr*650mm Proximal Release
	A:24, cup and ball:30, umbrella:50			

Remark: HC means half covered; FC means full covered.