


Spray Pipe

The solution for safe chromoendoscopy.

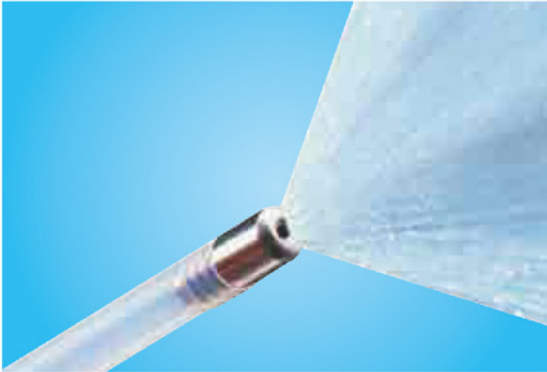

Intended use

- The staining of mucosa in the gastrointestinal tract contributes to an improved defining of details. Suspicious areas or tissue changes can be determined more reliable. The spray pipes help to applicate the staining solution precise.
- The innovative spray tip ensures efficient application of the liquid - either as an all-over spray by fine nozzles or as a concentrated jet for precise staining.

Special features

- High-quality manufacturing and kink resistant
- Easy to use and quick to deploy
- Smooth insertion of the catheter
- Efficient application of the solution through fine nozzles
- Two innovative spray heads available

Safe and easy insertion.

- The introduction of the spray pipe in the endoscope is easy and safe. Thanks to the kink-resistant material the catheter adapts at each turn of the endoscope. To use the staining solution the catheter is equipped with a Y-cannula.

Product code	Diameter of Jaws (mm)	Length (mm)	Units/Box
Straight shape			
CIT-PSG-Z-10x500	1.0	500	10
CIT-PSG-Z-10x1200	1.0	1200	10
CIT-PSG-Z-18x1200	1.8	1200	10
CIT-PSG-Z-18x1800	1.8	1800	10
CIT-PSG-Z-23x2300	2.3	2300	10
Mist shape			
CIT-PSG-W-10x500	1.0	500	10
CIT-PSG-W-10x1200	1.0	1200	10
CIT-PSG-W-18x1200	1.8	1200	10
CIT-PSG-W-18x1800	1.8	1800	10
CIT-PSG-W-23x2300	2.3	2300	10